

Brompheniramine

Waters Corporation




This is an Application Brief and does not contain a detailed Experimental section.

Abstract

This application brief demonstrates the analysis of brompheniramine using XTerra Columns.

Introduction

In this study the compound analyzed is brompheniramine.

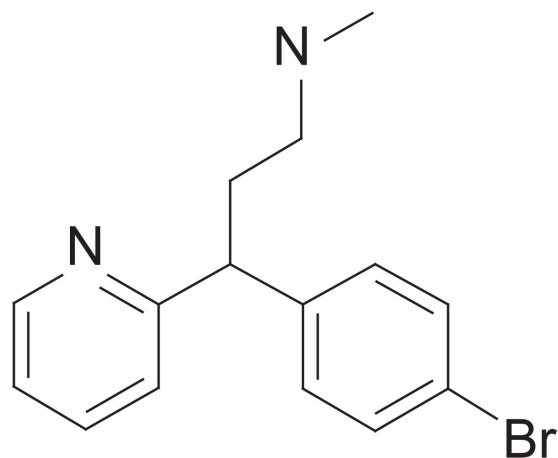

Brompheniramine


Experimental

Conditions

Column:	XTerra RP ₁₈ 4.6 x 150 mm, 5 μm
Part number:	186000492
Mobile phase:	pH 3.0: H ₂ O/ACN/100 mM NH ₄ COOH, pH 3.0 45:45:10 pH 7.0: H ₂ O/ACN/100 mM NH ₄ HCO ₃ , pH 7.0 50:40:10 pH 10.0 H ₂ O/ACN/100 mM NH ₄ HCO ₃ , pH 10.0 50:40:10

Flow rate:	1.0 mL/min
Injection volume:	5 μ L of 250 μ g/mL
Temperature:	30 $^{\circ}$ C
Detection:	UV @ 270 nm
Instrument:	Alliance 2695, 2996 PDA
Mobile Phase pH	USP Tailing
3.0	1.22
7.0	1.37
10.0	1.18

Results and Discussion

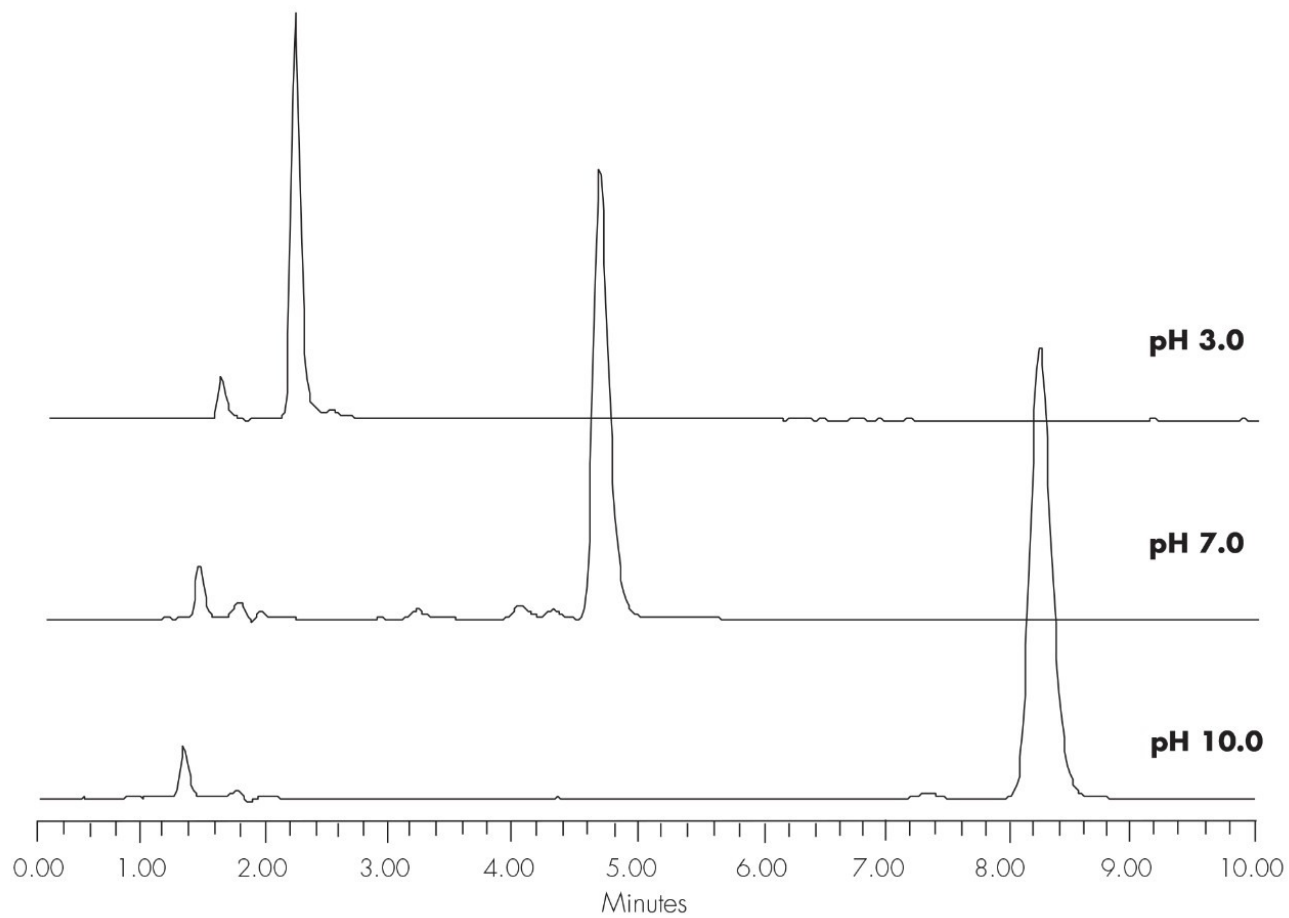

Featured Products

WA20738.018, June 2002



© 2021 Waters Corporation. All Rights Reserved.