

응용 자료

Pyrilamine

Waters Corporation




This is an Application Brief and does not contain a detailed Experimental section.

Abstract

This application brief highlights the analysis of pyrilamine using XTerra RP₁₈ columns.

Introduction

Pyrilamine has been analyzed in this application brief.

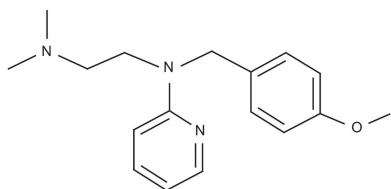

Pylamine

Experimental

HPLC Method

Column:	XTerra RP ₁₈ 4.6 x 150 mm, 5 μm (p/n: 186000492)
Mobile phase:	At pH 3.0: H ₂ O/ACN/100 mM NH ₄ COOH, pH 3.0, 45:45:10 At pH 7.0: H ₂ O/ACN/100 mM NH ₄ HCO ₃ , pH 7.0, 10:80:10 At pH 10.0: H ₂ O/ACN/100 mM NH ₄ HCO ₃ , pH 10.0, 50:40:10
Flow rate:	1.0 mL/min
Injection volume:	5 μL of 250 μg/mL
Temperature:	30 °C

Detection:


UV @ 280 nm

Instrument:

Alliance 2695, 2996 PDA

Mobile Phase pH	USP Tailing
3.0	1.28
7.0	1.23
10.0	1.09

Results and Discussion

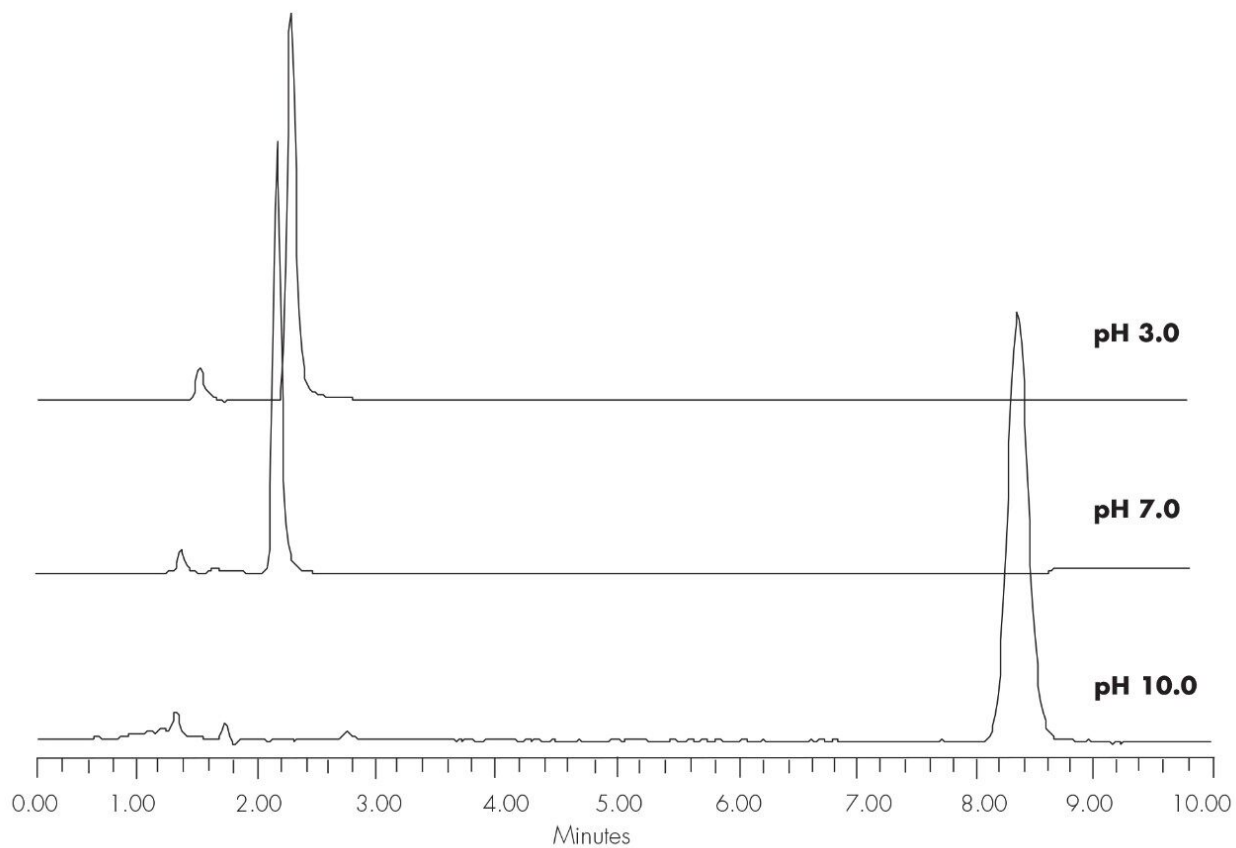

Featured Products

Alliance HPLC <<https://www.waters.com/514248>>

WA20738.094, June 2002

© 2021 Waters Corporation. All Rights Reserved.