

Dibucaine - LC/UV

Waters Corporation




This is an Application Brief and does not contain a detailed Experimental section.

Abstract

This application brief highlights the analysis of dibucaine using XTerra RP₁₈ columns.

Introduction

Dibucaine has been studied in this application brief.

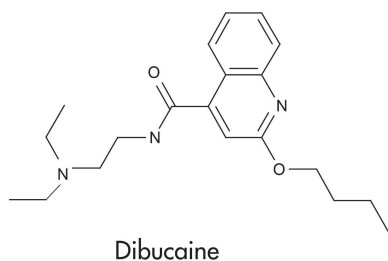

Experimental

HPLC Method

Column:	XTerra RP ₁₈ 4.6 x 150 mm, 5 μm (p/n: 186000492)
Mobile phase:	At pH 3.0: H ₂ O/ACN/100 Mm NH ₄ COOH, pH 3.0, 30:60:10 At pH 7.0: H ₂ O/ACN/100 mM NH ₄ HCO ₃ , pH 7.0, 35:55:10 At pH 10.0: H ₂ O/ACN/100 mM NH ₄ HCO ₃ , pH 10.0, 20:70:10
Flow rate:	1.0 mL/min
Injection volume:	5 μL of 250 μg/mL
Temperature:	30 °C

Detection:


UV @ 280 nm

Instrument:

Alliance 2695, 2996 PDA

Mobile phase pH	USP Trailing
3	1.23
7	1.19
10	1.18

Results and Discussion

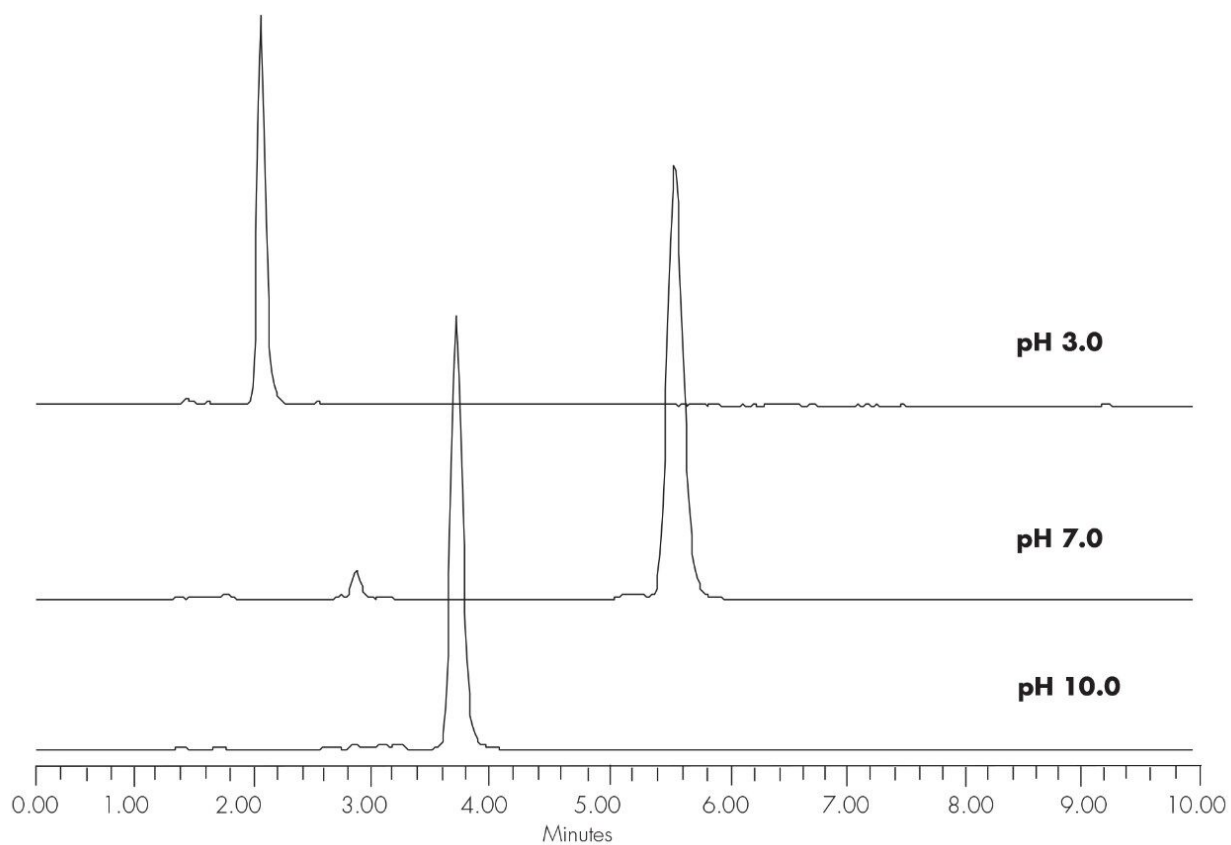

Featured Products

Alliance HPLC <<https://www.waters.com/514248>>

WA20738.034, June 2002

© 2021 Waters Corporation. All Rights Reserved.