


응용 자료

## Ranitidine (Symmetry C<sub>8</sub>)

---

Waters Corporation

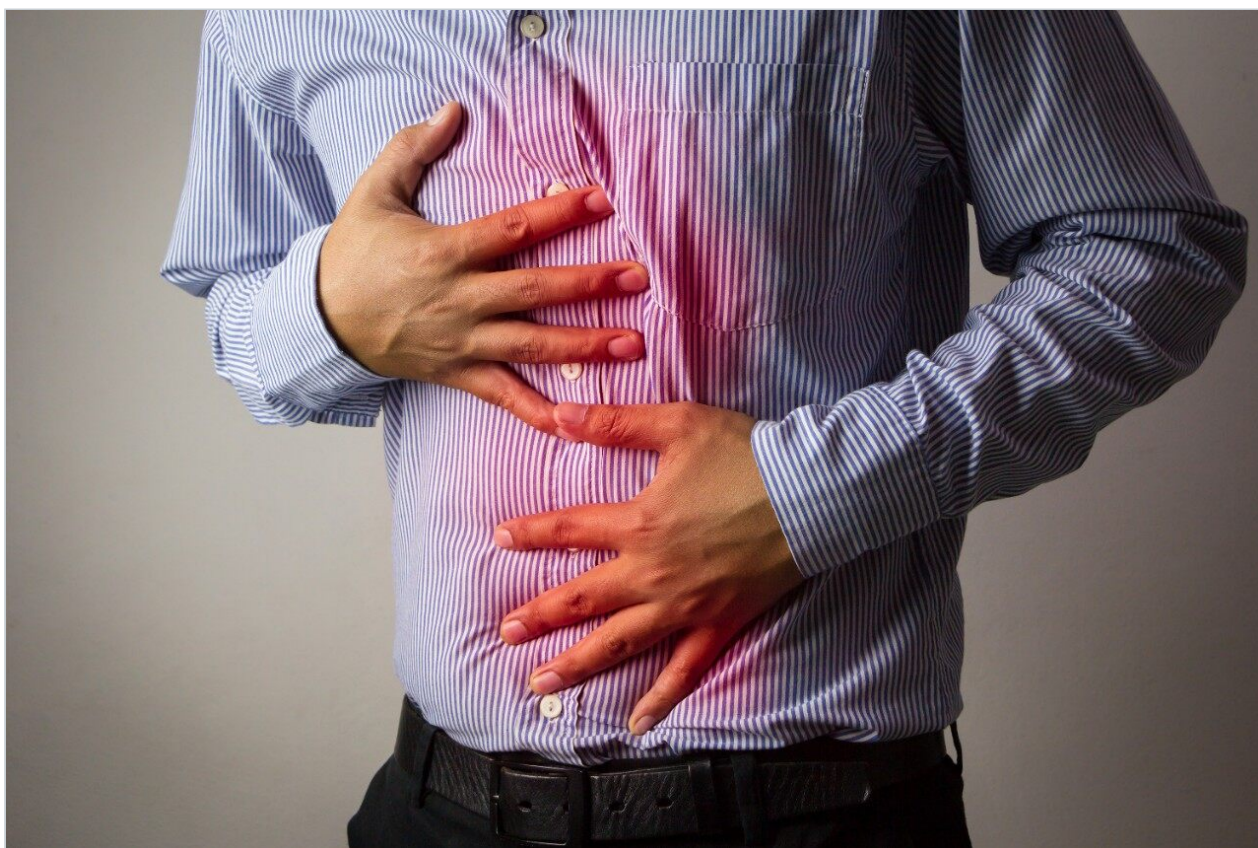

This is an Application Brief and does not contain a detailed Experimental section.

---


### Abstract

This application brief highlights the analysis of ranitidine using Symmetry C<sub>8</sub> Columns.

---

## Introduction

Ranitidine was analyzed using Symmetry C<sub>8</sub> Columns.

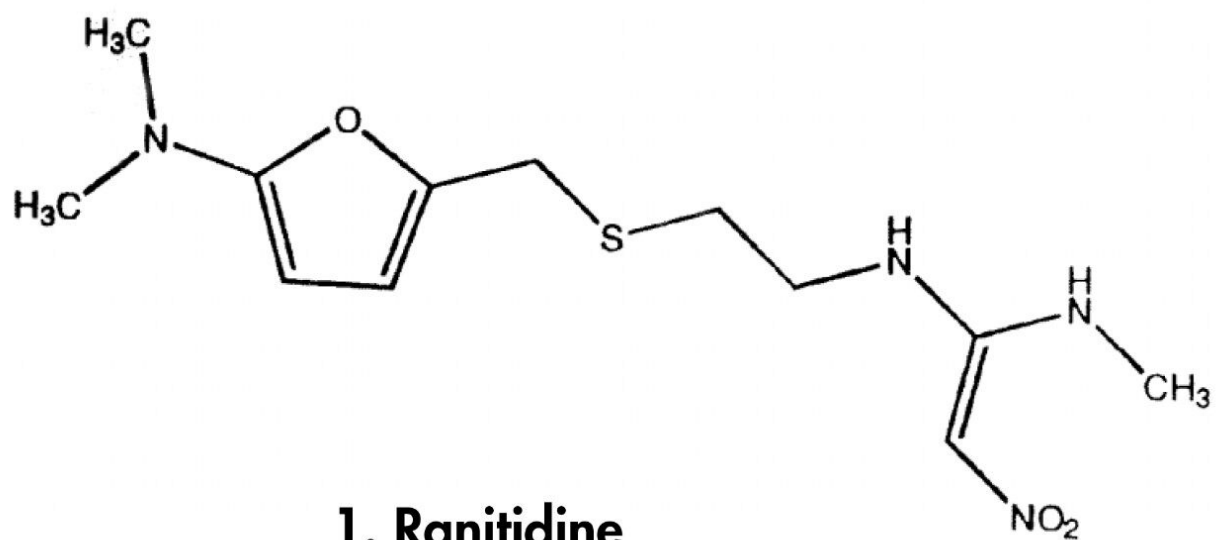

---

## Experimental

### HPLC Method

Column:	Symmetry C <sub>8</sub> , 3.9 x 150 mm, 5 μm
Guard column:	Symmetry Guard Column 3.9 x 20 mm, 5 μm
Part numbers:	Column - WAT046970, Guard - WAT054250
Mobile phase:	100 mM sodium acetate, 5mM camphosphonic acid, pH 6.2/ acetonitrile/methanol 80:10:10

Flow rate: 1.0 mL/min

Injection volume: 5  $\mu$ L of 100  $\mu$ g/ml sample


Detection: UV @ 254 nm

## USP Tailing Factor

1.112

---

## Results and Discussion

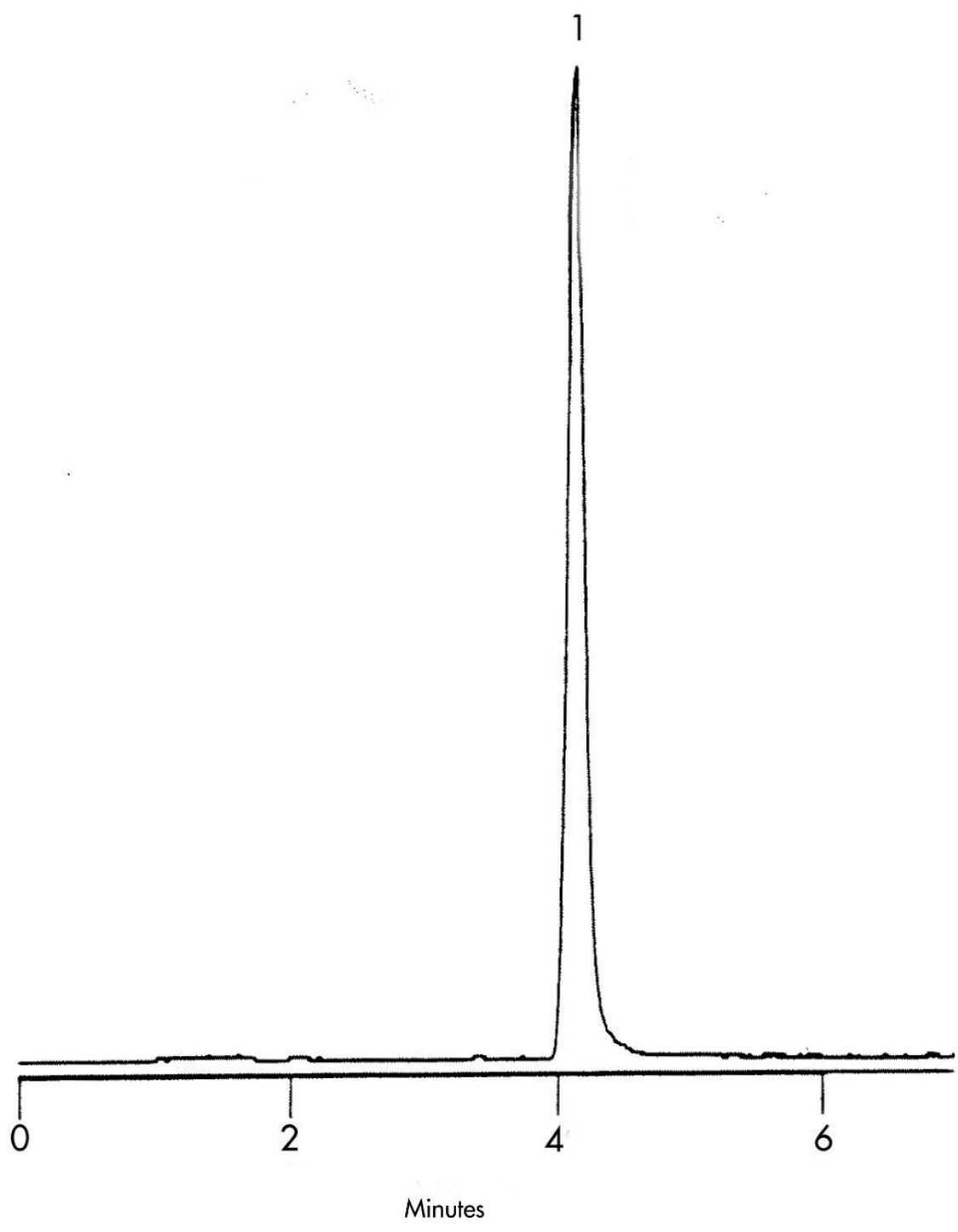

---

Featured Products

WA31763.143, June 2003

© 2021 Waters Corporation. All Rights Reserved.