

Note d'application

## Succinimides

---

Waters Corporation



This is an Application Brief and does not contain a detailed Experimental section.

---

### Abstract


This application brief demonstrates analysis of Succinimides.

---


## Introduction

The compounds used in this study are -


Compounds	USP Tailing Factors
Ethosuximide	0.96
Phensuximide	0.92
Methsuximide	0.91

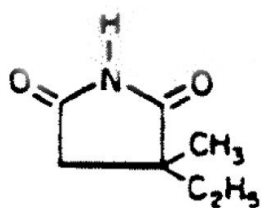

**1. Ethosuximide**

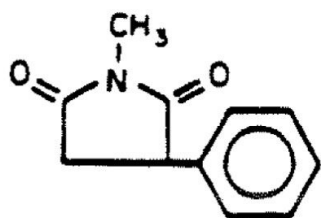

**2. Phensuximide**

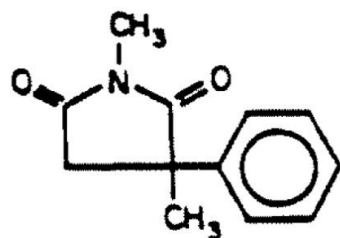

**3. Methsuximide**

---


## Experimental

### HPLC Method

Column:	Symmetry C <sub>8</sub> , 3.9 x 150 mm, 5 μm
Guard column:	Sentry Guard Column 3.9 x 20 mm, 5 μm
Part numbers:	Column - WAT046970, Guard - WAT054250
Mobile phase:	100 mM potassium phosphate, pH 6.9/methanol/water 20:30:50
Flow rate:	1.0 mL/min
Injection volume:	10 μL of 40–500 μg/mL sample
Detection:	UV @ 214 nm

---

## Results and Discussion

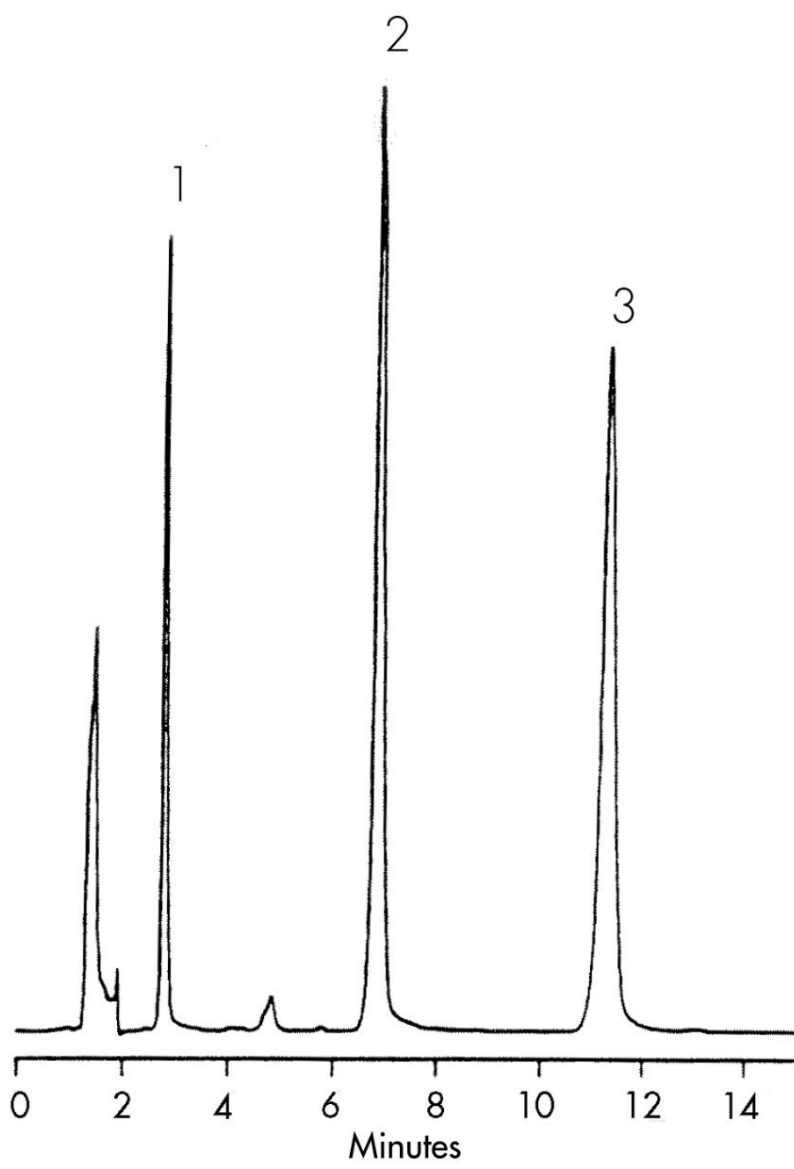

---

## Featured Products

WA31763.153, June 2003

© 2021 Waters Corporation. All Rights Reserved.