


Trazodone

Waters Corporation

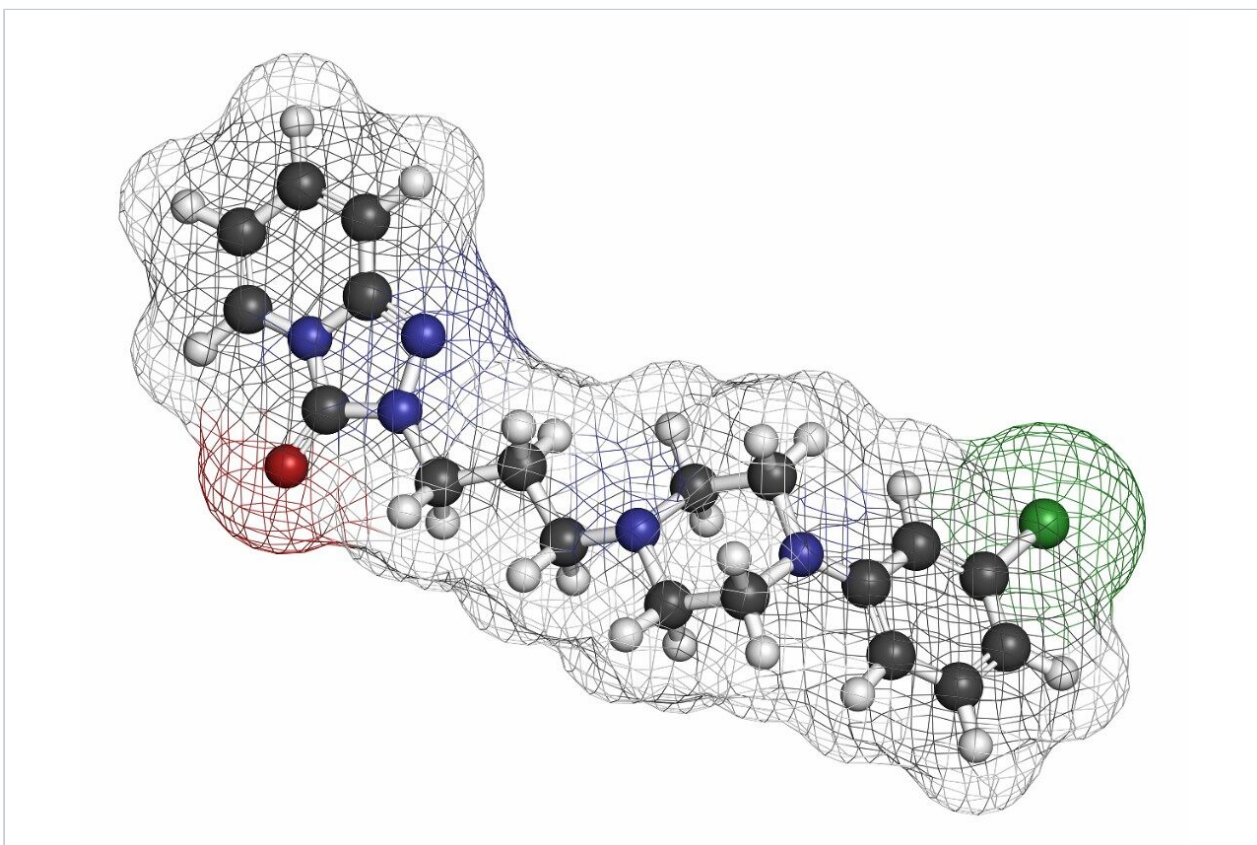

This is an Application Brief and does not contain a detailed Experimental section.


Abstract

This application brief demonstrates analysis of trazodone.


Introduction

The compounds used in this study are –


Compounds	USP Tailing Factors
1-(3-chlorophenyl)peperazine	1.31
Maprotiline (I.S.)	1.21
Trazodone	1.12

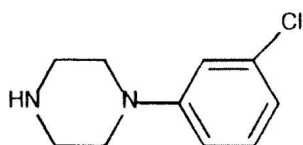

1. 1-(3-chlorophenyl) peperazine

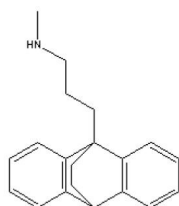

2. Maprotiline (I.S.)

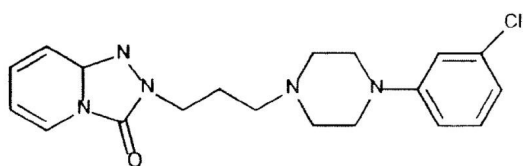


3. Trazodone

Experimental

HPLC Method

Column:	Symmetry Shield RP ₈ , 3.9 x 150 mm, 5 μm
Part numbers:	WAT200655
Mobile phase:	50 mM phosphate, pH7.0/methanol 36:64
Flow rate:	1.0 mL/min
Injection volume:	25 μL
Detection:	UV @ 250 nm

Results and Discussion

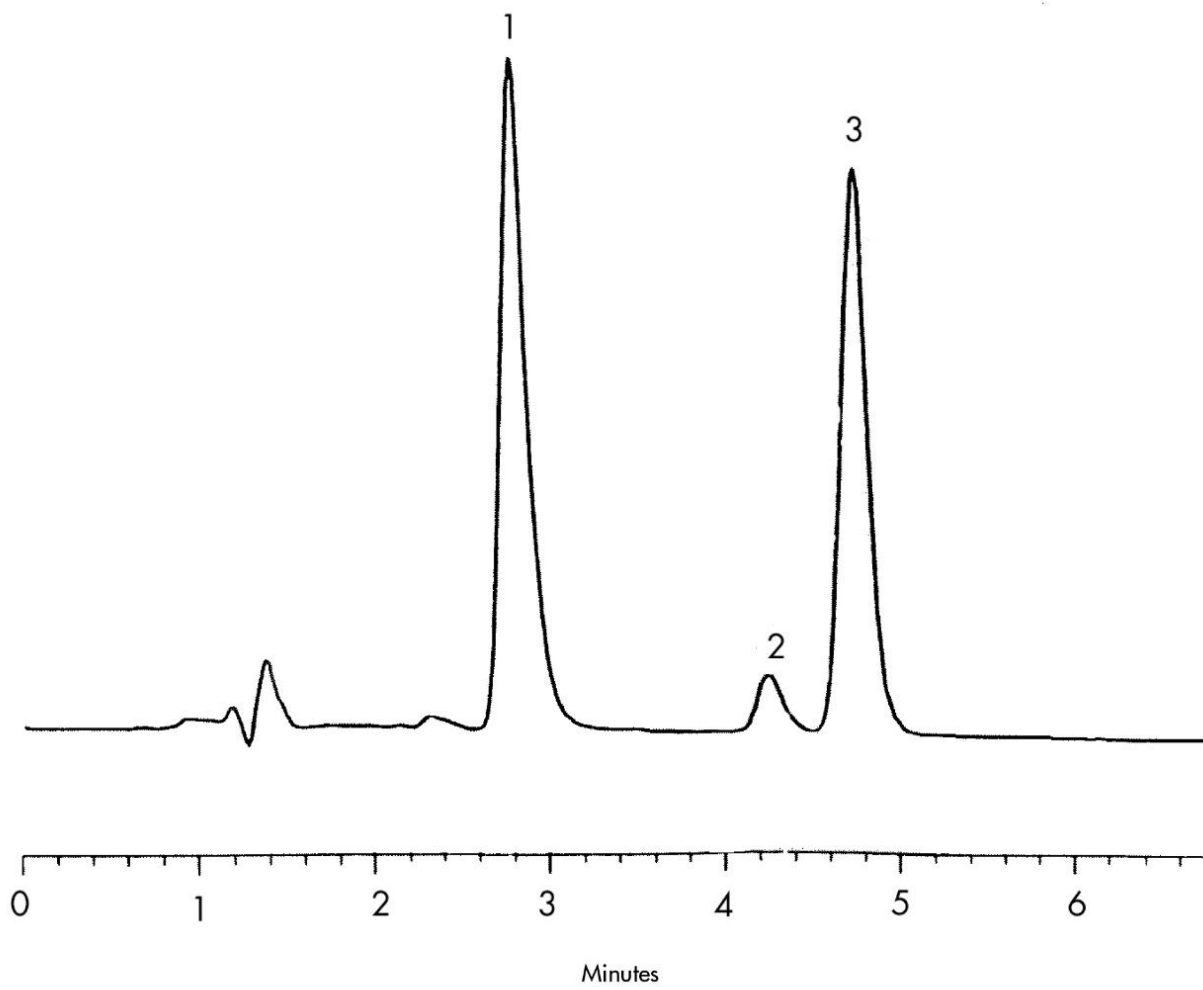

Featured Products

WA31763.171, June 2003

© 2021 Waters Corporation. All Rights Reserved.