

Application Note

## Betamethasone

---

Waters Corporation



This is an Application Brief and does not contain a detailed Experimental section.

---

### Abstract


This application brief demonstrates analysis of betamethasone.

---


## Introduction

The compounds used in this study are –

Compounds	USP Tailing Factors
Betamethasone	1.12
Propylparaben (I.S.)	1.09

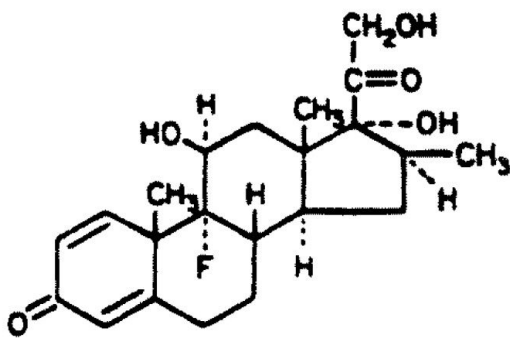

**1. Betamethasone**

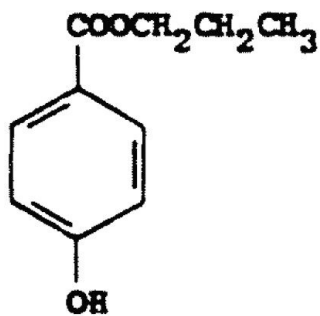

**2. Propylparaben  
(internal standard)**

---


## Experimental

### HPLC Method

Column:	Symmetry C <sub>8</sub> , 3.9 x 150 mm, 5 μm
Guard column:	Sentry Guard Column 3.9 x 20 mm, 5 μm
Part numbers:	Column - WAT046970, Guard - WAT054250
Mobile phase:	Water/acetonitrile/phosphoric acid 63:37:0.1
Flow rate:	1.0 mL/min
Injection volume:	10 μL of 50 μg/mL sample
Detection:	UV @ 254 nm

---

## Results and Discussion

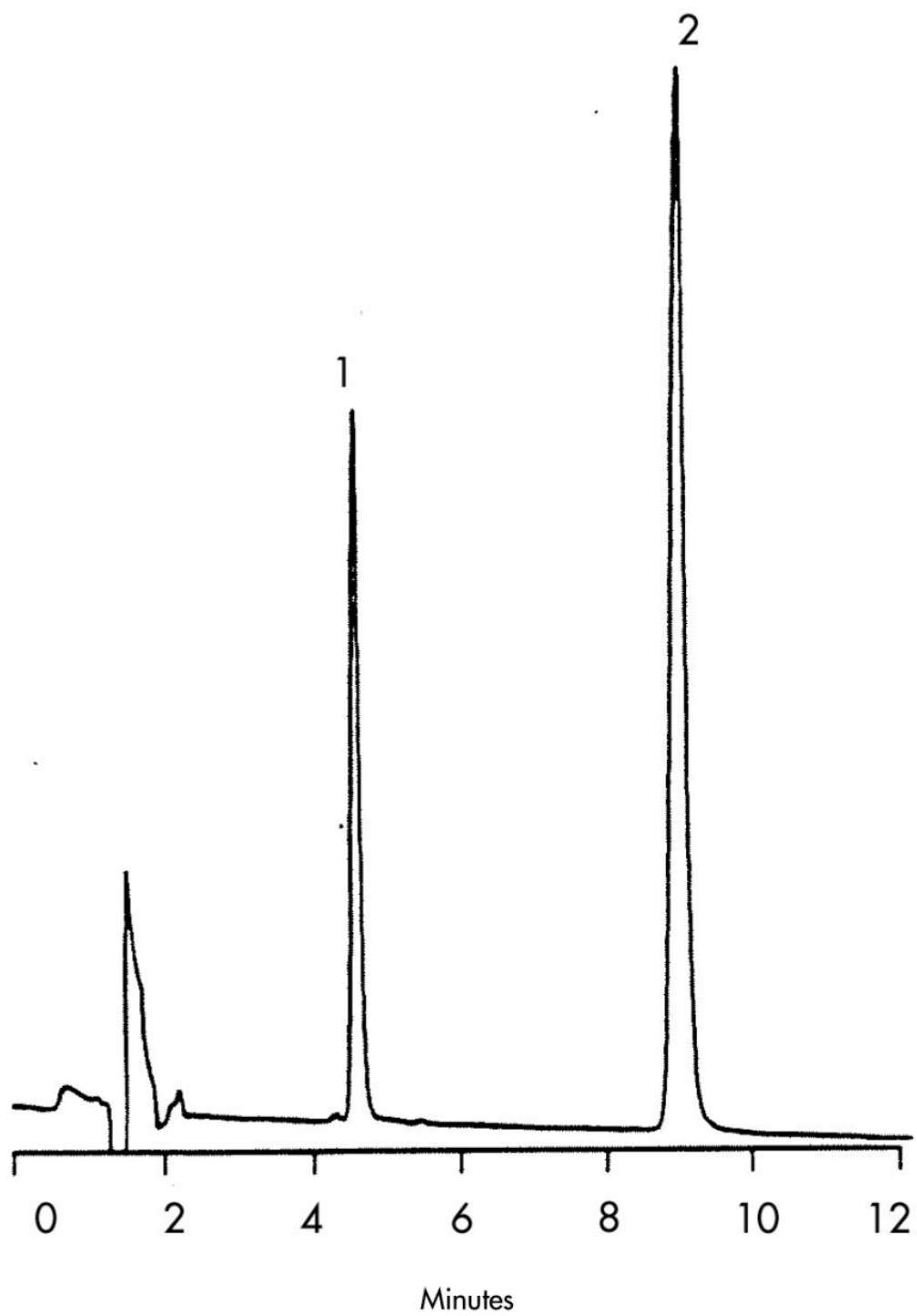

---

Featured Products

WA31763.37, June 2003

© 2021 Waters Corporation. All Rights Reserved.