

Application Note

## Betamethasone in Syrup Formulation

---

Waters Corporation



This is an Application Brief and does not contain a detailed Experimental section.

---

### Abstract


This application brief demonstrates analysis of betamethasone in syrup formulation.

---


## Introduction

The compounds used in this study are -

Compounds	USP Tailing Factors
Betamethasone	1.14
Propylparaben (I.S.)	1.12

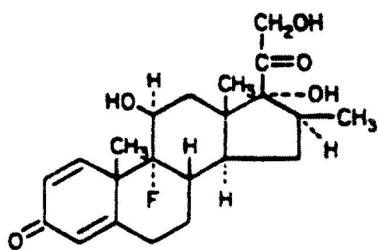

1. Betamethasone

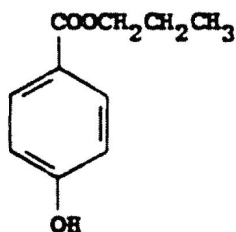

2. Propylparaben (I.S.)

---


## Experimental

### HPLC Method

Column:	Symmetry C <sub>8</sub> , 3.9 x 150 mm, 5 μm
Guard column:	Sentry Guard Column 3.9 x 20 mm, 5 μm
Part numbers:	Column - WAT046970, Guard - WAT054250
Mobile phase:	Water/acetonitrile/phosphoric acid 63:37:0:1
Flow rate:	1.0 mL/min
Injection volume:	10 μL of 150 μg/mL sample
Detection:	UV @ 254 nm

---

## Results and Discussion

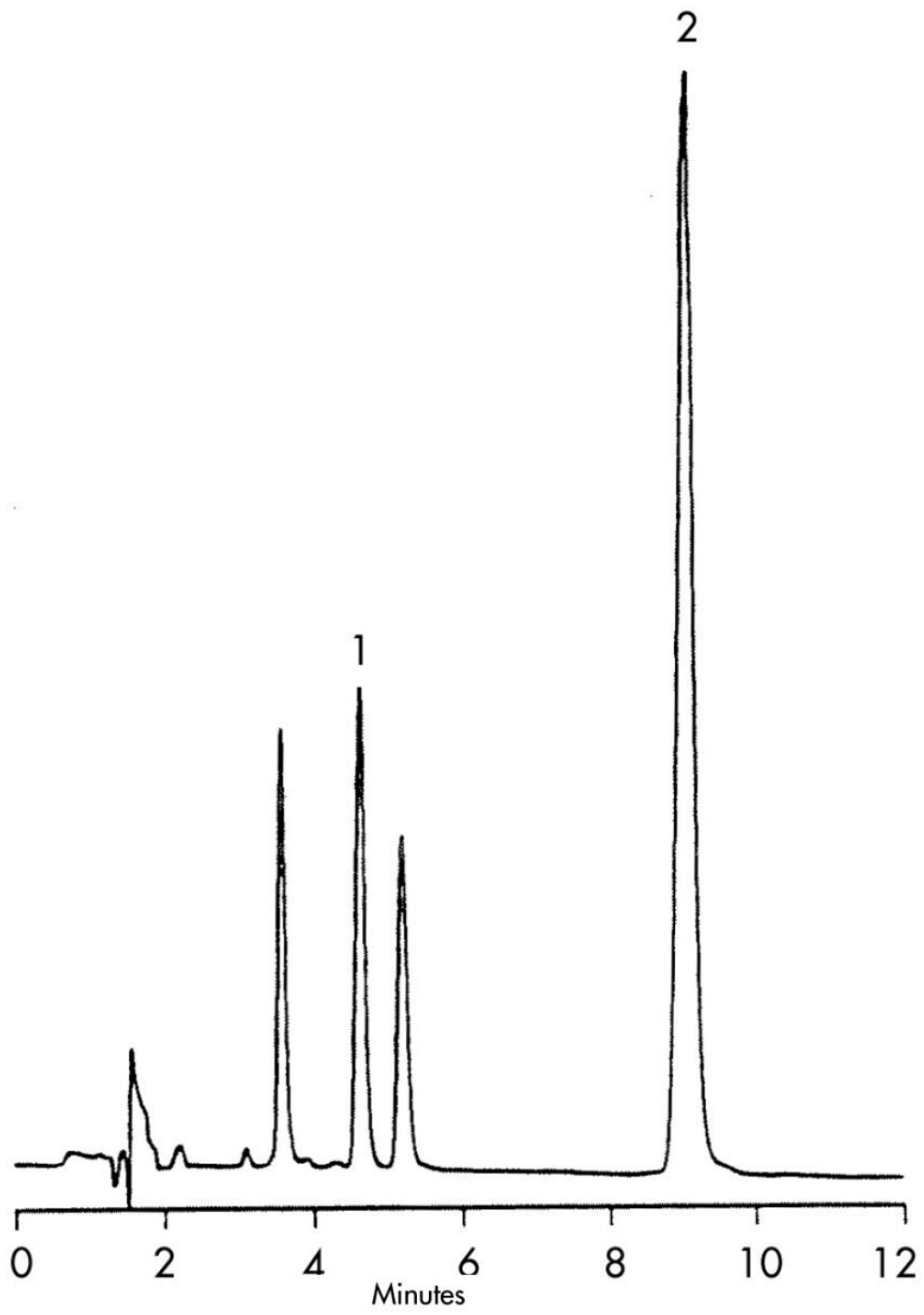

---

Featured Products

WA31763.38, June 2003

© 2021 Waters Corporation. All Rights Reserved.