


Nota de aplicación

Carbadox

Waters Corporation

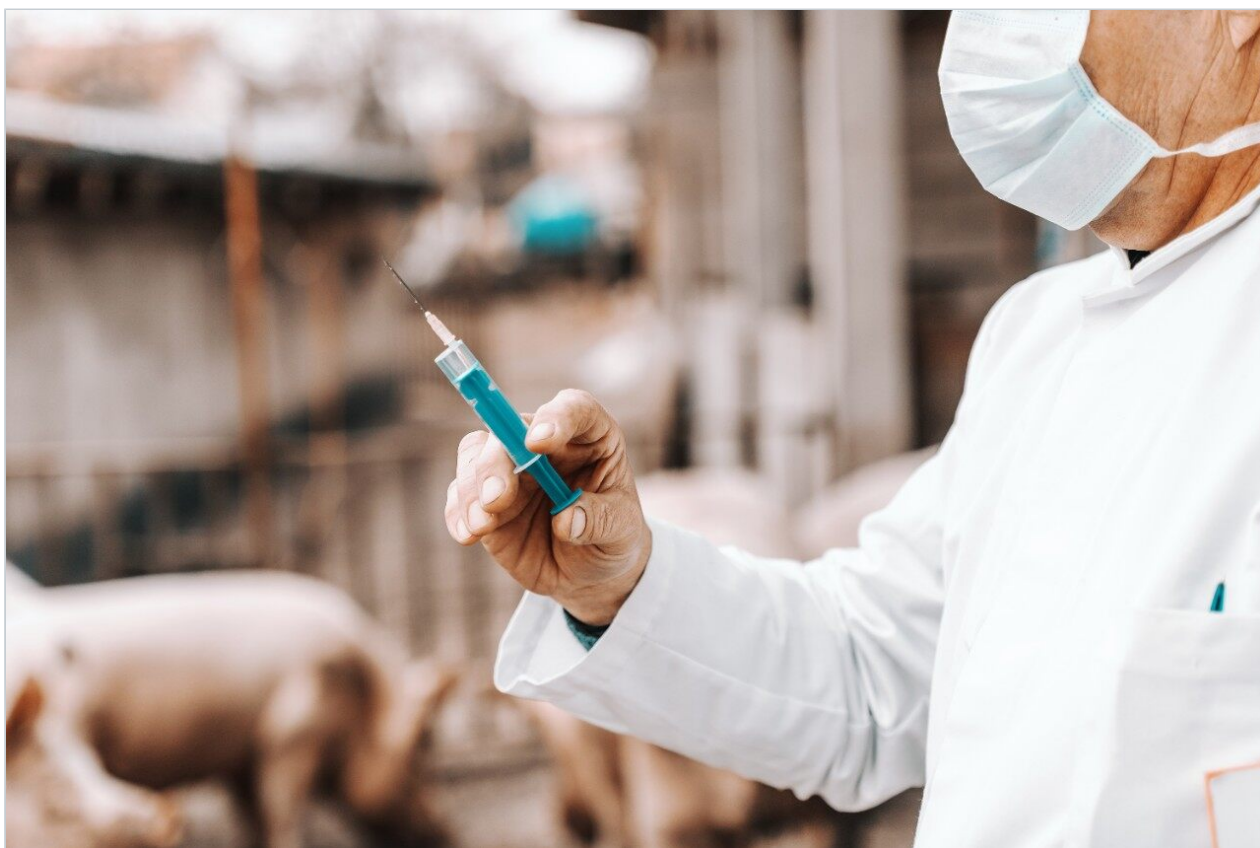


This is an Application Brief and does not contain a detailed Experimental section.

Abstract

This application brief demonstrates analysis of carbadox.

Introduction

The compound analyzed in this study is carbadox (USP Tailing Factor - 1.1).

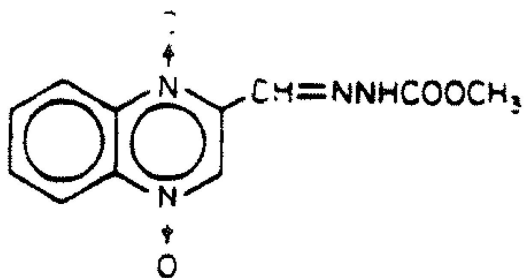

1. Carbadox

Experimental


HPLC Method

Column:	Symmetry C ₁₈ , 4.6 x 100 mm, 3.5 μm
Part number:	WAT066220
Mobile phase:	20 mM ammonium acetate, pH 5.0/Acetonitrile 85:15
Flow rate:	1.4 mL/min
Injection volume:	15 μL of 0.025 mg/mL
Temperature:	30 °C

Detection:

UV @ 254 nm

Results and Discussion

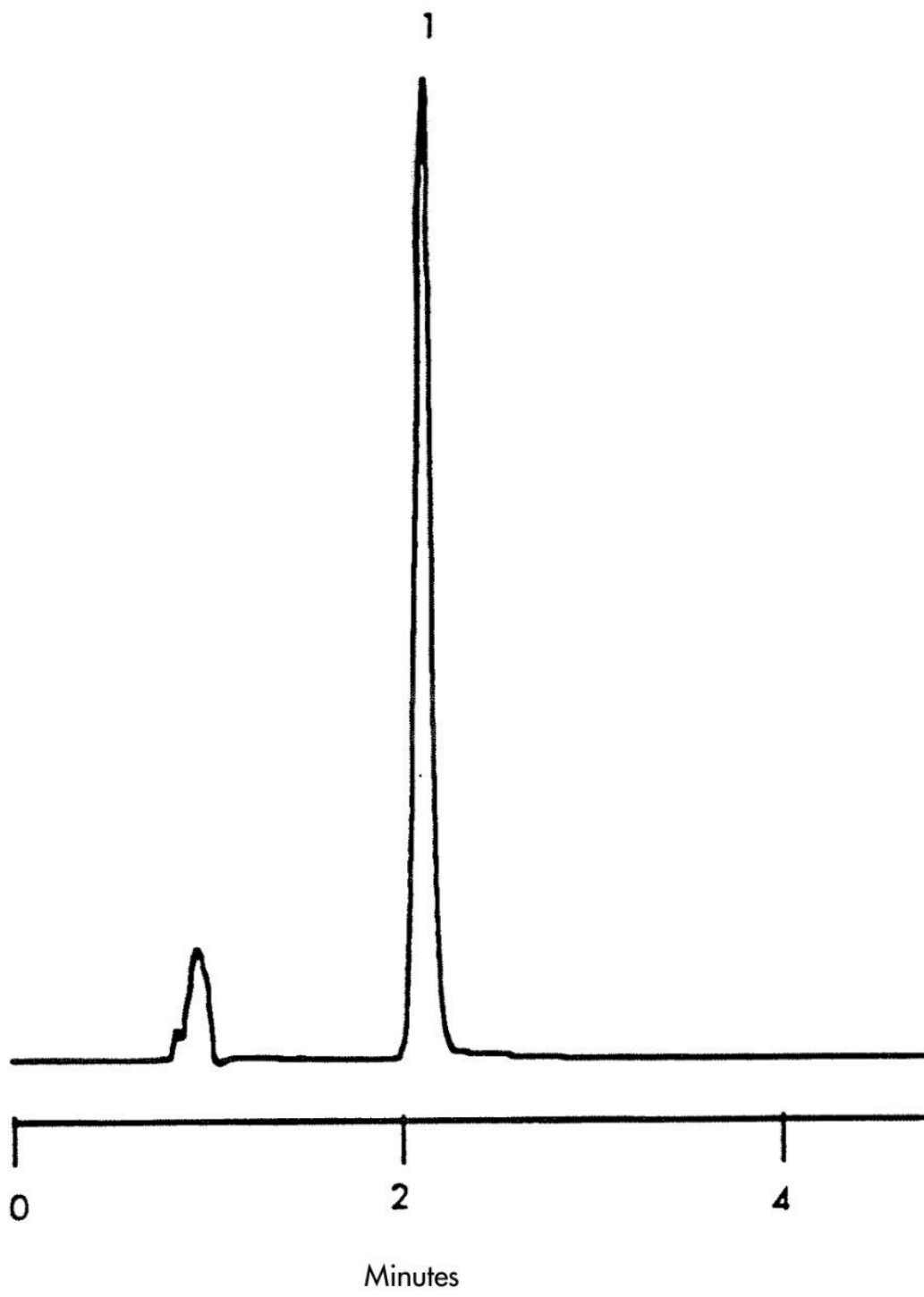

Featured Products

WA31763.47, June 2003

© 2021 Waters Corporation. All Rights Reserved.