

应用纪要

Diltiazem HCL

Waters Corporation


This is an Application Brief and does not contain a detailed Experimental section.

Abstract

This application brief demonstrates the analysis of diltiazem HCL using Symmetry Columns.

Introduction

The compound analyzed in this study is diltiazem HCL.

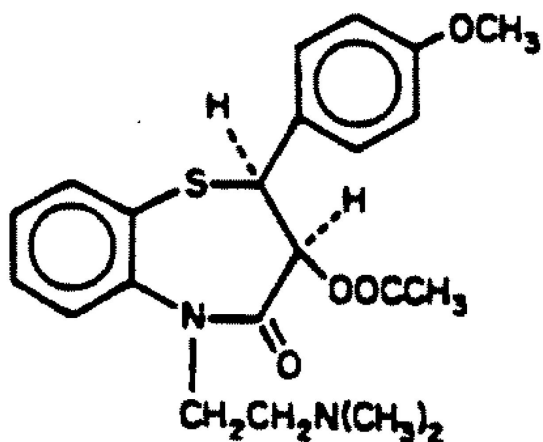

1. Diliazem HCL

Experimental

HPLC Method

| | |
|-------------------|---|
| Column: | Symmetry C ₈ , 3.9 x 150 mm, 5 μm |
| Guard Column: | Symmetry Guard Column 3.9 x 20 mm, 5 μm |
| Part numbers: | Column - WAT046970, Guard - WAT054250 |
| Mobile phase: | 100 mM sodium acetate, 5mM camphosulfonic acid, pH 6.2/acetonitrile/methanol 40:30:30 |
| Flow rate: | 1.0 mL/min |
| Injection volume: | 10 μL of 50 μg/mL sample |


Detection:

UV @ 254 nm

USP Tailing Factors

1. 1.08

Results and Discussion

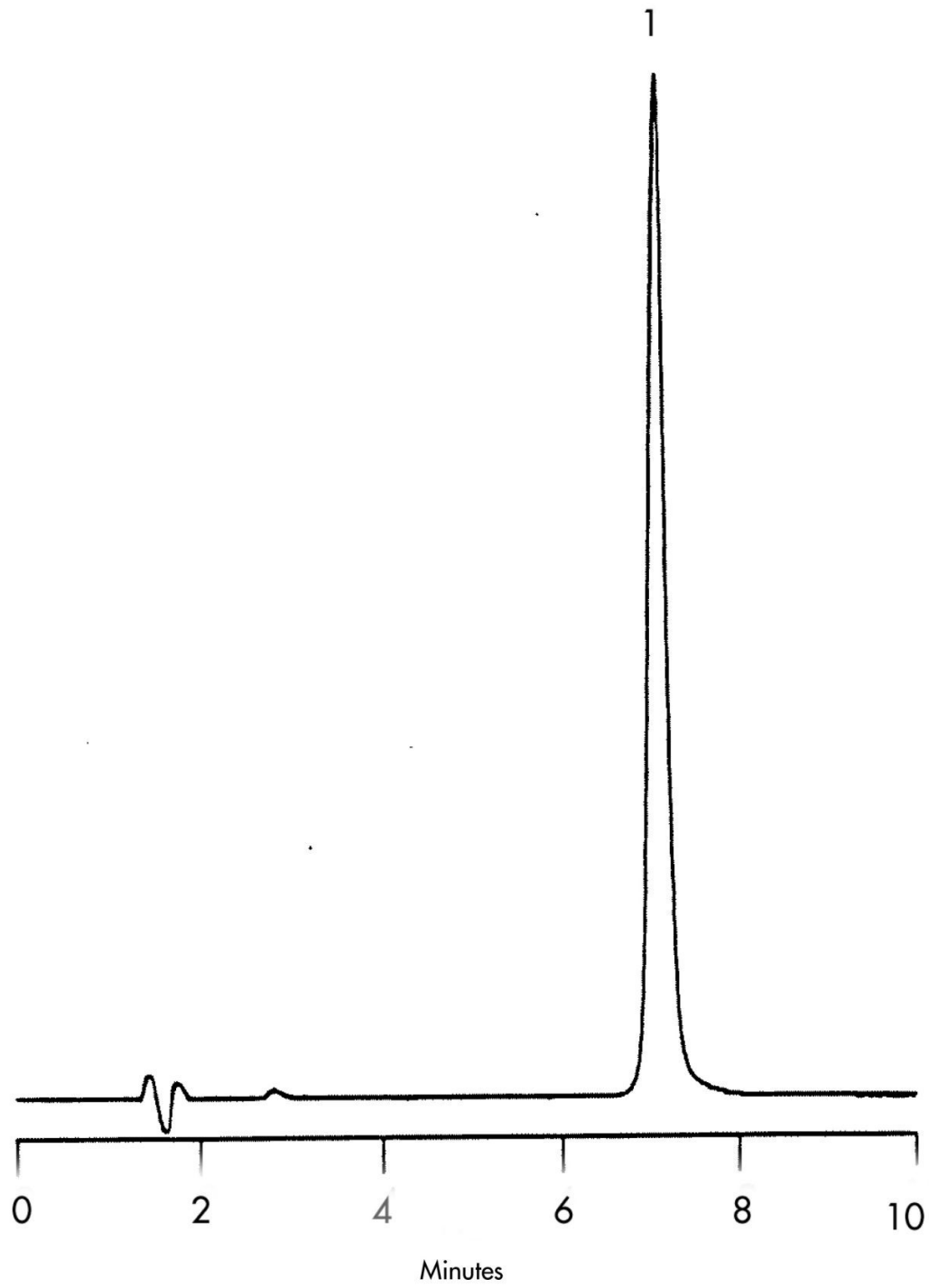

Featured Products

WA31763.75, June 2003

© 2022 Waters Corporation. All Rights Reserved.