

## Separation of a Tryptic Digest of Phosphorylase B

---

Waters Corporation



This is an Application Brief and does not contain a detailed Experimental section.

---

### Abstract

This application brief demonstrates separation of a tryptic digest of phosphorylase B.

---

## Introduction

The compound used in this study is tryptic digest of phosphorylase B.


---

## Experimental

Column:	BioSuite PA Series (2.1 x 150 mm)
Eluent A:	0.02% TFA in water
Eluent B:	0.016% TFA in acetonitrile
Gradient:	0–56% B in 70 minutes
Flow rate:	0.16 mL/min
Injection volume:	15 $\mu$ L
Temperature:	40 °C
Detection:	UV 220 nm

---

## Results and Discussion

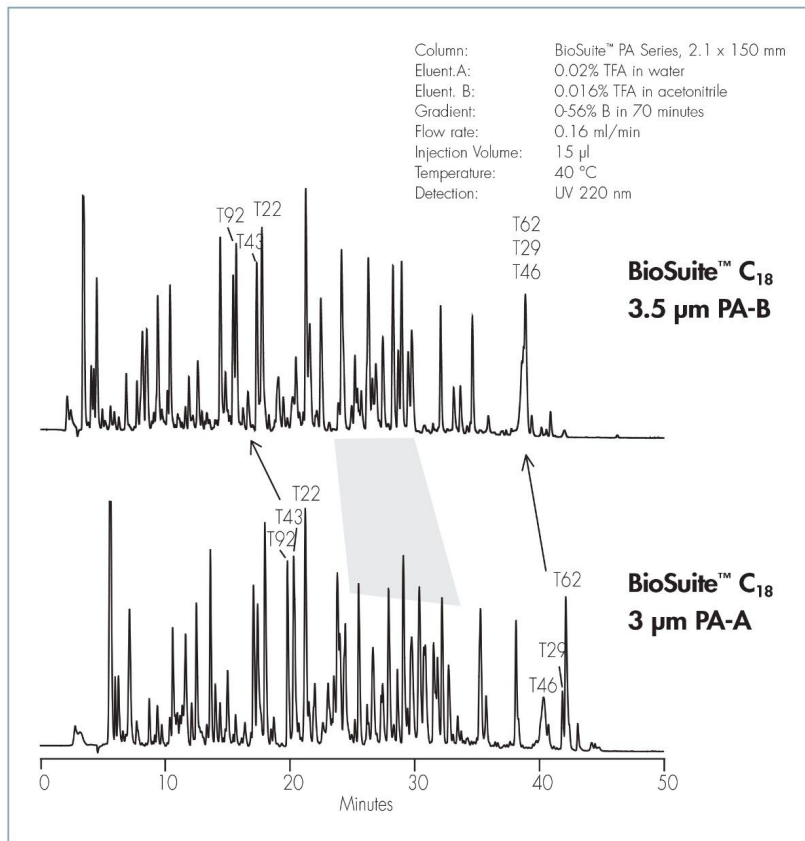

## Featured Products

WA31792, June 2003